



# **DataView v4 User Guide**

Software version 4

---

## DataView v4 Software User Guide

The contents of this DataView™ User Guide are the property of Respiratory Technology Corporation (Restech). No part of this publication may be reproduced, stored in a retrieval system, transmitted by any means – electronic, mechanical, photocopy, recording or otherwise be used in any form, by any means other than its intended purpose, without the prior written permission of Respiratory Technology Corporation.

Respiratory Technology Corporation  
114 Holmes Rd. Suite 201  
Houston, TX 77045  
www.restech.com  
858-673-3700 Ph 858-673-3783 Fx

Toll Free Customer Service and Technical Support  
U.S. Customers: 800-352-1512  
International Customers: +1 (858) 673-3700

This publication is supplied for use with the Dx–pH Measurement System™ only. Refer to the latest version of the Dx–pH Measurement System Instructions for Use in order to properly use the DataView v4 software with the Dx–pH Measurement System.

Respiratory Technology Corporation is responsible for the reliability and performance of the Dx–pH Measurement System™ only if the following conditions are met:

- All repairs and modifications are carried out by Respiratory Technology Corporation or its assigned agents.
- The equipment is used in accordance with the Instructions for Use.

**CAUTION: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.**

The Dx-pH Measurement System complies with the requirements of Health Canada Medical Devices Regulation SOR/98-282, in accordance with ISO 13485:2003 Medical Devices – Quality Management Systems – Requirements for Regulatory Purposes, under the Canadian Medical Devices Conformity Assessment System (CMDCAS), and TGA Australian regulatory guidelines for medical devices (ARGMD).

CE 0843 – The device complies with the Medical Device Directive (MDD, 93/42/EEC). The four-digit number refers to the Notified Body who issued the certificate relative to this directive.



EMERGO EUROPE  
Prinsessegracht 20  
2514 AP The Hague  
The Netherlands

# DataView v4 Software User Guide

## Table of Contents

<b>1. Introduction to DataView v4</b>	<b>45</b>
<b>2. System Requirements</b>	<b>6</b>
<b>3. Starting the DataView Software</b>	<b>7</b>
<b>4. Main Navigation Menu Selection “New”</b>	<b>8</b>
4.1 Entering and Saving New Study Information	8
4.1.1 Patient Details Tab	9
4.1.2 Custom Symptoms Tab	11
4.1.3 Study Details Tab	12
<b>5. Main Navigation Menu Selection “Retrieve”</b>	<b>14</b>
5.1 Retrieve a Completed Study	14
<b>6. Main Navigation Menu Selection “Open”</b>	<b>15</b>
6.1 Open an Existing Study	15
<b>7. Main Navigation Menu Selection “Print”</b>	<b>18</b>
7.1 Print Graph	18
7.2 Print Report	20
<b>8. Main Navigation Menu Selection “Graph”</b>	<b>22</b>
8.1 View and Edit pH Study Graph	22
8.2 View and Edit Event Graph	25
8.3 Graph Toolbar	26
8.3.1 Add Event and Symptom Information	27
8.3.2 Remove Event and Symptom Information	28
8.3.3 Delete Event Information	29
8.3.4 Undo Changes	29
8.3.5 Zoom In	30
8.3.6 Zoom Out	30
8.3.7 Y-Axis (pH Scale) Adjustment	31
8.3.8 Delete Meals	33

# DataView v4 Software User Guide

<b>9. Main Navigation Menu Selection “Report”</b>	<b>34</b>
9.1 View and Edit Reports	34
9.1.1 Site Information Block	35
9.1.2 Restech Information Block	35
9.1.3 Patient Study Information Block	36
9.1.4 pH Study Graph	36
9.1.5 Acid Exposure Graph	37
9.1.6 Event Calculations Table	38
9.1.7 Symptom Correlation Table	39
9.1.8 RYAN Score Table	39
9.1.9 Notes	40
9.1.10 Signature	40
<b>10. Main Navigation Menu Selection “Save”</b>	<b>41</b>
<b>11. Main Navigation Menu Selection “Settings”</b>	<b>42</b>
11.1 Settings	42
11.1.1 General Tab	42
11.1.2 Normal Acid Exposure	43
11.1.3 Colors Tab	44
11.1.4 Event Settings	45
11.1.5 Path to Data	45
11.2 Help	45
<b>12. Uninstalling DataView v4 Software</b>	<b>46</b>

# DataView v4 Software User Guide



## 1. Introduction to DataView v4

DataView is provided by Restech authorized personnel on media (thumb drive, SD card or disk) and installed in person during a field support visit; or during an online technical support session.

DataView features the following enhancements:

- 1) STANDARD CLOCK  
Displays time in 24 hour clock format.
- 2) SIMPLIFIED MENU AND TABS  
Main Navigation Menu and dialog boxes for user-friendly control.
- 3) EDITABLE PH STUDY GRAPH  
Allows the user to precisely edit the pH Study Graph using the Graph Toolbar.
- 4) COMPREHENSIVE REPORT  
Displays a comprehensive Report page with the pH Study Graph, Event Calculations Table, Symptom Correlations Table, Acid Exposure Graph, and RYAN Score.
- 5) MULTIPLE BASELINES  
Displays event and symptom calculations for four different baseline pH levels.
- 6) FAST EVENT AND SYMPTOM CALCULATIONS  
Automatically calculates and presents events and symptoms in an easily viewable table on the Report page.
- 7) ADJUSTABLE COLORS  
Allows the user to customize graph color settings.

DataView has a simplified Graphical User Interface that incorporates an easily accessible menu.

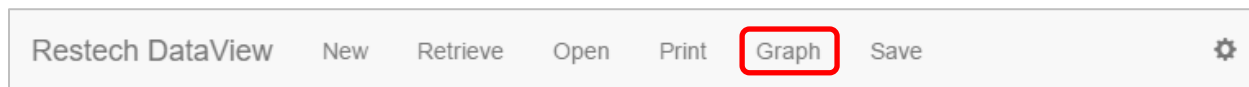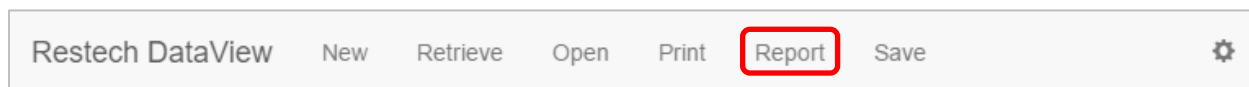

Report and Graph toggle allows the user to easily switch between Report and Graph view using the Main Navigation Menu.

## 2. System Requirements

### Operating System Compatibility

The DataView program runs on the following operating systems: Windows 7®, Windows 8®, and Windows 10®.

### Minimum System Requirements

The following are the minimum recommended system requirements for computers using DataView:

Description	Minimum Requirement
Processor type	Intel Pentium® Dual-Core or equivalent
Random Access Memory	512 MB
Free Disk Space	150 MB for program; 6 MB per study file
Pointing device	Mouse, touchpad, or equivalent device
Data Card reader	Compatible with SD memory card

### Program/User Security Requirements

Entry into the DataView software program is not password protected. Users are expected to password protect the computing device on which the program runs. Access to patient data and files will be controlled by users.

### User permissions


The user must be logged into an account to properly install DataView.

### 3. Starting the DataView Software

To start DataView, click on the Restech DataView Icon from the Taskbar or Desktop.



DataView launches with the Main Navigation Menu and Settings in the DataView Workspace.

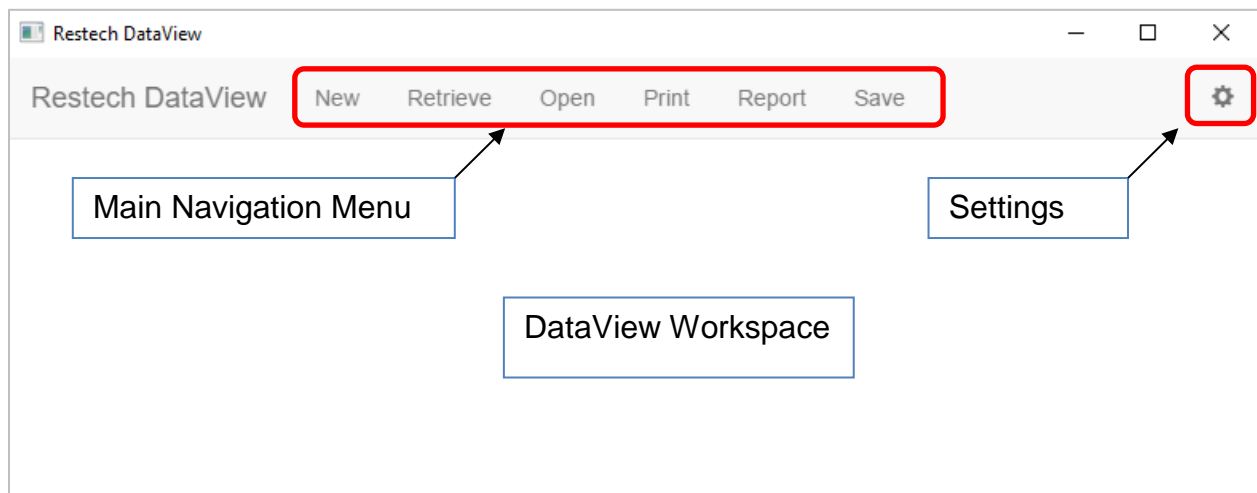


## 4. Main Navigation Menu Selection “New”


### 4.1 Entering and Saving New Study Information

To set-up a New Study, insert a RESTECH formatted SD card into your computer's SD card reader. (If a navigation pane emerges with the data card file directory, close it and continue.)

The "New" Main Navigation Menu Selection displays a New Study dialog box with input fields for the user to enter Patient Details, Custom Symptoms, and Study Details in separate tabs.

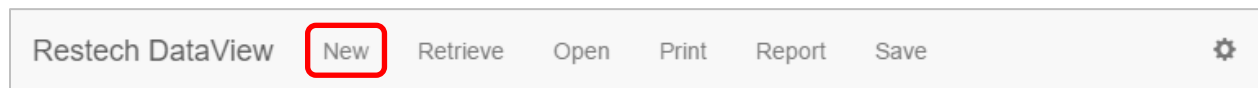

Click on "New" and enter Patient Details, Custom Symptoms, and Study Details according to the following **Sections 4.1.1 – 4.1.3**.

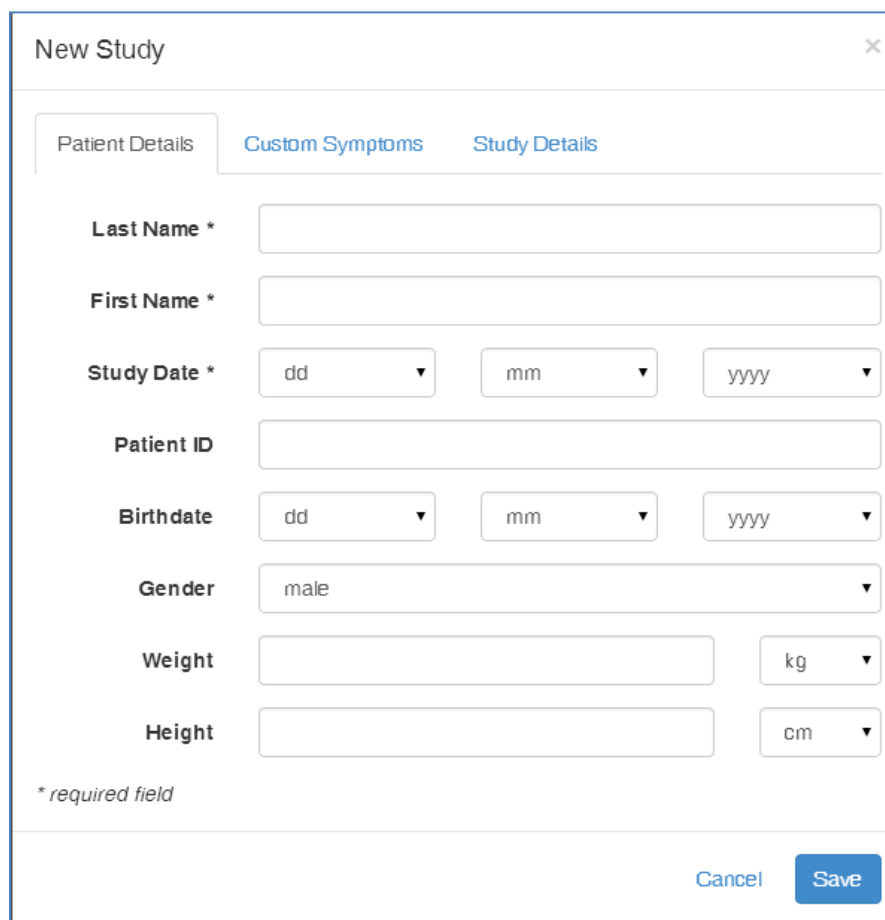

New Study

Patient Details Custom Symptoms Study Details

Last Name \*

First Name \*

Study Date \*

Patient ID

Birthdate

Gender

Weight

Height

\* required field

Cancel Save



## DataView v4 Software User Guide

### 4.1.1 Patient Details Tab

Click on the "New" Main Navigation Menu Selection to open a New Study dialog box. There are 3 tabs provided for entering information. The first tab displays Patient Details input fields.

Click in the field boxes to input patient data according to the available controls:

Field	Available Controls
Patient Last Name	30 characters
Patient First Name	30 characters
Study Date	dd-mm-yyyy [selectable, drop down]
Patient ID	30 characters
Birth Date	dd-mm-yyyy [selectable, drop down]
Gender	male/female [selectable, toggle]
Patient Weight	n pounds (n kg) [selectable, toggle]
Patient Height	n inches (n cm) [selectable, toggle]

New Study

Patient Details

Custom Symptoms

Study Details

Last Name \*

First Name \*

Study Date \*

dd

mm

yyyy

Patient ID

Birthdate

dd

mm

yyyy

Gender

male

Weight

kg

Height

cm

\* required field

Cancel

Save

## DataView v4 Software User Guide

The required fields Last Name, First Name, and Study Date are labeled with an asterisk. If the user clicks "Save" without entering patient information into these field boxes, the required fields turn red.

The screenshot shows a 'New Study' form with three tabs: 'Patient Details', 'Custom Symptoms', and 'Study Details'. The 'Patient Details' tab is active. A red callout box labeled 'Required Information' points to the 'Last Name \*', 'First Name \*', and 'Study Date \*' fields. The 'Study Date' field consists of three dropdown menus for 'dd', 'mm', and 'yyyy'. Other fields include 'Patient ID', 'Birthdate' (with 'dd', 'mm', 'yyyy' dropdowns), 'Gender' (with a 'male' dropdown), 'Weight' (with a 'kg' dropdown), and 'Height' (with a 'cm' dropdown). A note at the bottom left states '\* required field'. At the bottom right are 'Cancel' and 'Save' buttons.

**Required Information**

**Last Name \***

**First Name \***

**Study Date \***

**Patient ID**

**Birthdate**

**Gender**

**Weight**

**Height**

\* required field

Cancel Save

Enter patient information into the required fields and any other desired fields.

Information entered in Patient Details will appear on the Report header after the study is completed.

Next, click the Custom Symptoms tab to enter specific symptom information.

## DataView v4 Software User Guide

### 4.1.2 Custom Symptoms Tab

Click on the "New" Main Navigation Menu Selection to open a New Study dialog box. The second tab displays Custom Symptoms input fields.

Click in the Symptom fields to input symptoms according to the available controls:

Field	Available Controls
Symptom 1	4 characters
Symptom 2	4 characters
Symptom 3	4 characters

The screenshot shows a 'New Study' dialog box with three tabs: 'Patient Details', 'Custom Symptoms', and 'Study Details'. The 'Custom Symptoms' tab is active, displaying three input fields labeled 'Symptom 1', 'Symptom 2', and 'Symptom 3'. Each field is a text box with a light gray border. At the bottom right of the dialog, there are two buttons: 'Cancel' (light blue) and 'Save' (dark blue). A close button (X) is located in the top right corner of the dialog box.

The symptoms will appear on the Recorder screen during the pH study.

Next, click the Study Details tab to enter study specific information.

## DataView v4 Software User Guide

### 4.1.3 Study Details Tab

Click on the "New" Main Navigation Menu Selection to open a New Study dialog box. The third tab displays Study Details input fields. Enter information according to the available controls:

Field	Available Controls
Technician	30 characters
Physician	30 characters
Referring	30 characters
Study Type	pharyngeal [selectable]
Indication	Infinite characters
Antireflux Med.	30 characters

The screenshot shows a 'New Study' dialog box with three tabs: 'Patient Details', 'Custom Symptoms', and 'Study Details'. The 'Study Details' tab is active, showing input fields for 'Technician', 'Physician', 'Ref. Physician', 'Study Type' (a dropdown menu with 'pharyngeal' selected), 'Indication', and 'Antireflux Med.'. At the bottom right, there are 'Cancel' and 'Save' buttons.


Information entered in Study Details will appear on the Report header after the study is completed.

See the next section for Saving the New Study information.

## DataView v4 Software User Guide

### 4.1.4 Save New Study Information


Click "Save" from any tab in the New Study dialog box. The 'Please wait Saving' box is displayed.

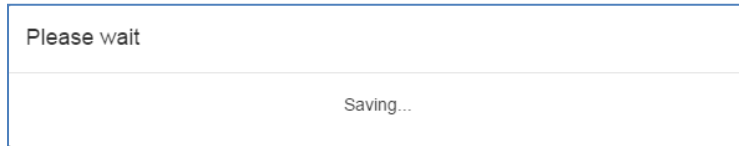

Once the New Study information is saved to the RESTECH SD card, the 'Please wait Saving' box disappears and the DataView Workspace is displayed.

Note: You can click "Save" on any of the three tabs and still enter additional information on the other tabs. Just click "Save" again to save the additional information.


If there is an existing study on the RESTECH SD card, when the user clicks "Save" from any tab in the New Study dialog box, a 'File exists Overwrite?' dialog box is displayed.

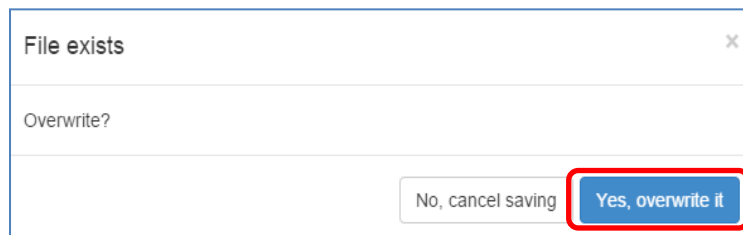

If the existing study has already been saved to the computer, click "Yes, overwrite it." The 'Please wait Saving' box is displayed.

Once the New Study information is saved to the RESTECH SD card, the 'Please wait Saving' box disappears and the DataView Workspace is displayed.


If the existing study has not been saved to the computer, click "No, cancel saving." The 'File exists Overwrite?' dialog box disappears and the DataView Workspace is displayed. To save the existing study on the SD card follow the instructions in Section 5.1 Retrieve a Completed Study, then return to Section 4.1 Entering and Saving New Study Information.

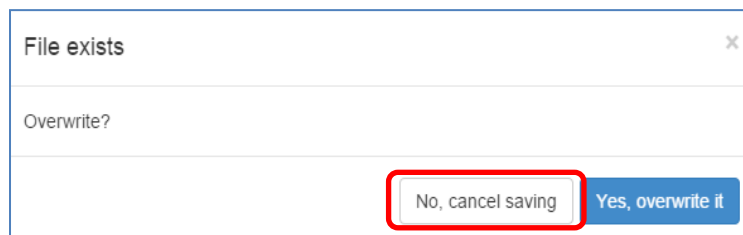

### 5. Main Navigation Menu Selection “Retrieve”

#### 5.1 Retrieve a Completed Study


The "Retrieve" Main Navigation Menu Selection allows the user to save a completed study into the Restech data folder and to pull up the study data into the Graph Workspace.

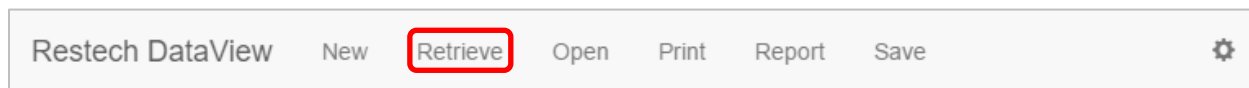

To retrieve a patient study, remove the RESTECH formatted SD card from the Dx-Recorder and insert it into your computer's SD card reader. You may use the SD card reader supplied by Restech if your computer does not have a built-in card reader. A navigation pane may display the data card file directory. Close it and continue.


Start the DataView software by clicking on the Restech DataView Icon. DataView launches with the Main Navigation Menu in the blank DataView Workspace.

Click on the "Retrieve" Main Navigation Menu Selection. A 'Please wait Retrieving' dialog box is displayed.

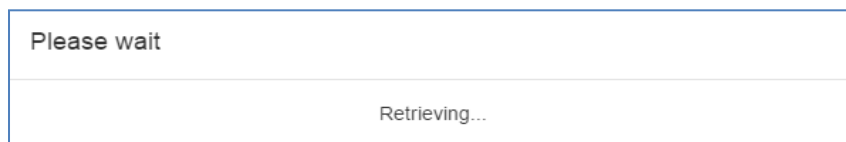

When DataView completes the data retrieval, the 'Please wait Retrieving' dialog box will disappear and the Graph Workspace will show the patient's pH Study Graph.


If there is no patient study on the SD card, a 'Sorry No PatientData.txt file was found' dialog box is displayed. Click "OK" to return to the DataView Workspace. Insert another SD card with a patient study to retrieve a study.

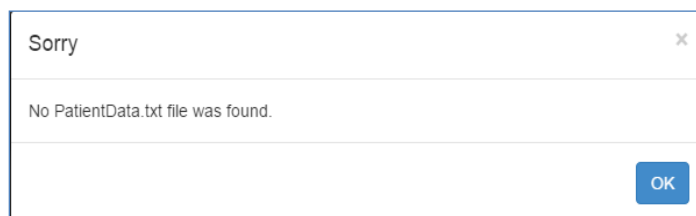

## 6. Main Navigation Menu Selection “Open”


### 6.1 Open an Existing Study

The "Open" Main Navigation Menu Selection displays an Open Study dialog box with a file search field. A list of recently saved patient files is displayed. A "Browse" button is provided for the user to open a patient study from anywhere on the computer or a shared drive location.

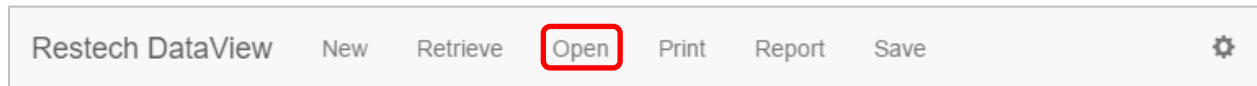

To open a patient study, click on the "Open" Main Navigation Menu Selection.


With the initial installation of the DataView software, there will be two example patient study files available in the Open Study list, 'Heartburn Joey' and 'Really Normal'. Additional patient study files will be added when completed studies have been retrieved according to Section 5.1.

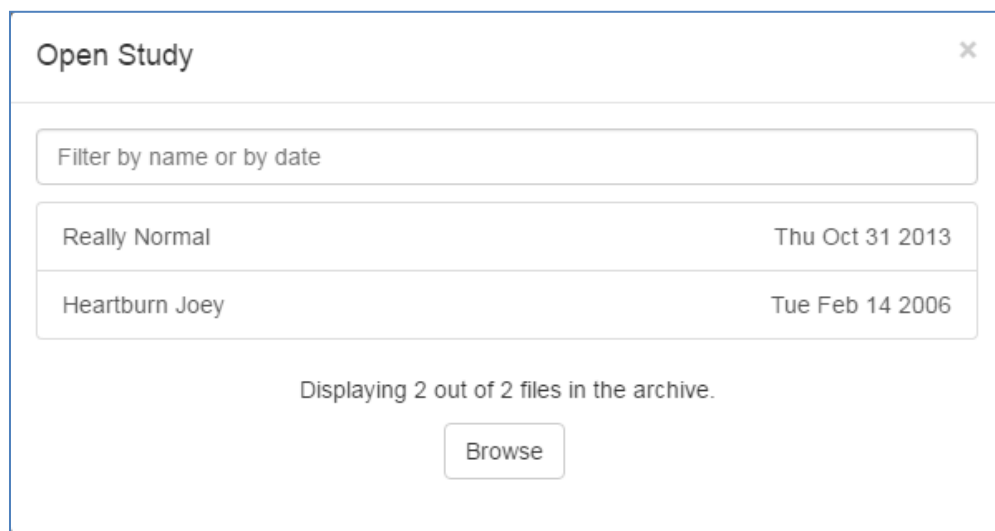

Once 10 or more additional patient study files have been retrieved and saved in DataView, the Open Study dialog box will display a list of the 10 most recent patient study files in the format 'Last name First name study date.' The dialog box includes a searchable list by name or date.

## DataView v4 Software User Guide


Click on a name or study date to open the file. DataView will open the pH study Graph.

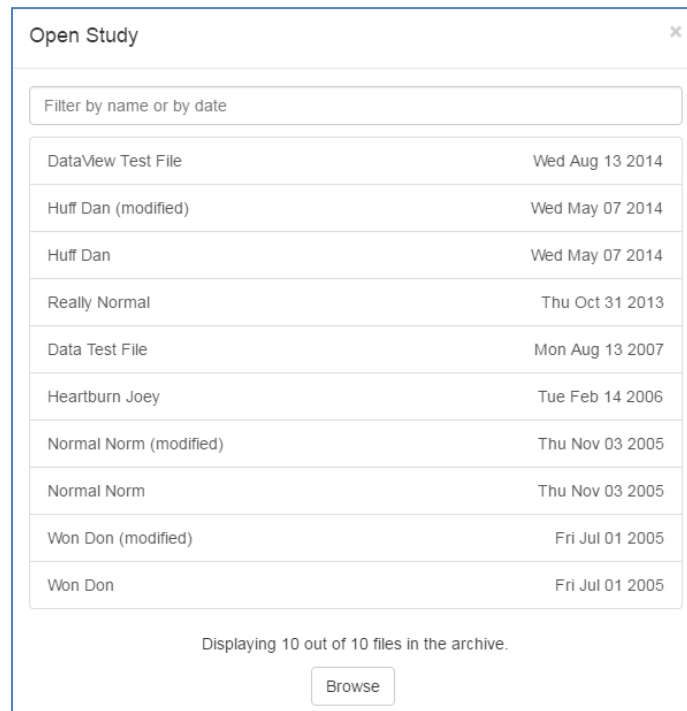

The 'Open Study' dialog box contains a search field labeled 'Filter by name or by date'. Below it is a table with 10 rows, each containing a file name and a date. At the bottom, it says 'Displaying 10 out of 10 files in the archive.' and has a 'Browse' button.


File Name	Date
DataView Test File	Wed Aug 13 2014
Huff Dan (modified)	Wed May 07 2014
Huff Dan	Wed May 07 2014
Really Normal	Thu Oct 31 2013
Data Test File	Mon Aug 13 2007
Heartburn Joey	Tue Feb 14 2006
Normal Norm (modified)	Thu Nov 03 2005
Normal Norm	Thu Nov 03 2005
Won Don (modified)	Fri Jul 01 2005
Won Don	Fri Jul 01 2005

If the patient's name or study date does not appear in the list, enter the first few letters of the patient name or the study date in the search field "Filter by name or by date". When the desired file is displayed, click on the file and DataView will open the pH study Graph.

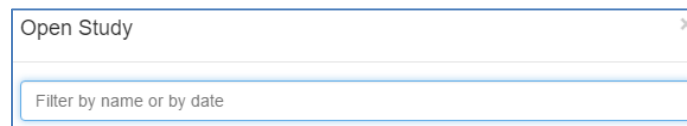

The 'Open Study' dialog box is shown with the search field 'Filter by name or by date' highlighted by a blue border.

To open a patient study that is not in the Open Study list and is not available in the search field, click the "Browse" button.

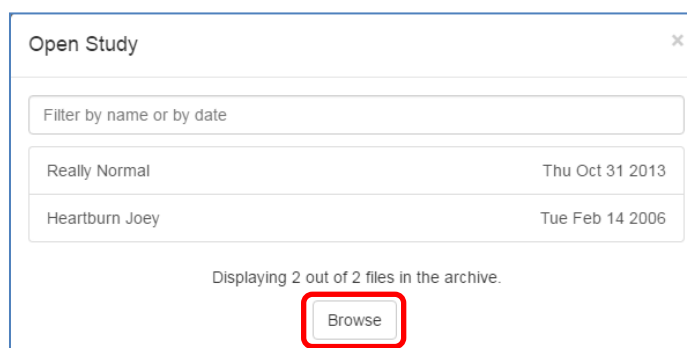


The 'Open Study' dialog box is shown with the 'Browse' button highlighted by a red rectangle. The search field contains the text 'Filter by name or by date'. The table below shows two files: 'Really Normal' (Thu Oct 31 2013) and 'Heartburn Joey' (Tue Feb 14 2006). At the bottom, it says 'Displaying 2 out of 2 files in the archive.'

File Name	Date
Really Normal	Thu Oct 31 2013
Heartburn Joey	Tue Feb 14 2006




## DataView v4 Software User Guide

A Windows Navigation Pane will be displayed. This allows navigation to the computer or a shared drive location to find and open the desired patient study. Select the desired file name and click Open.

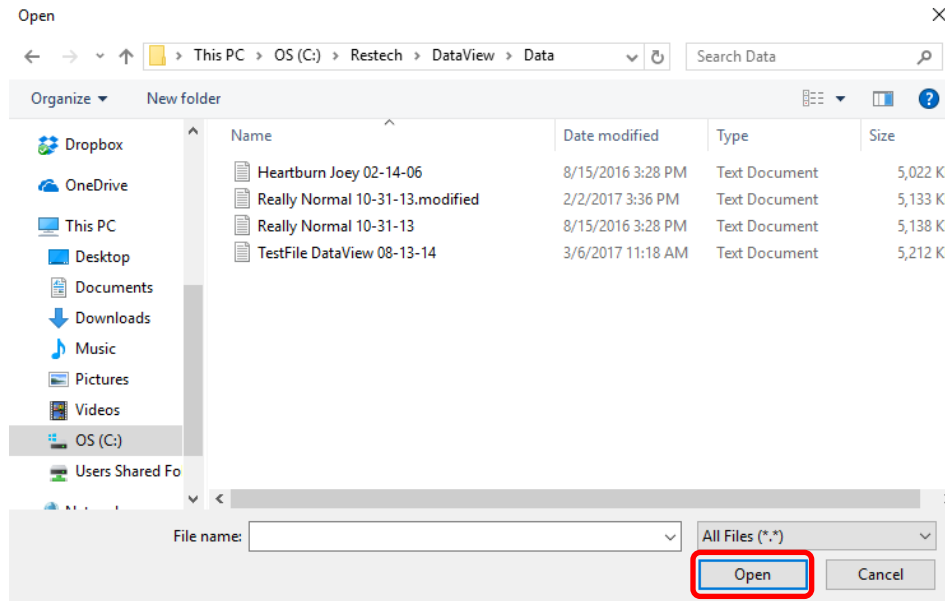

DataView will display "Loading..." at the bottom of the Open Study dialog box.

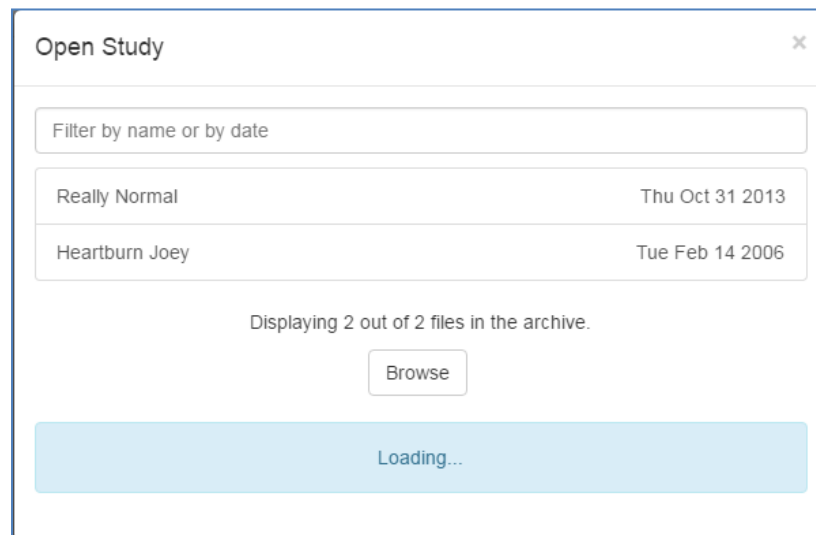


The Open Study dialog box will disappear and DataView will open the pH study Graph.

### 7. Main Navigation Menu Selection “Print”


#### 7.1 Print Graph

The "Print" Main Navigation Menu Selection allows the user to save a PDF document of Graphs and Reports.

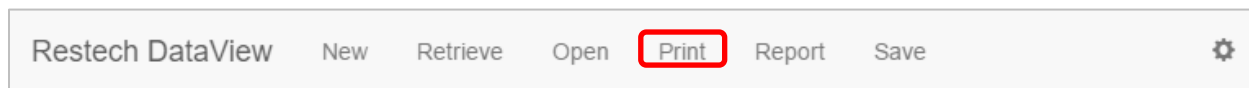

To print a Graph from the Graph Workspace, select "Print" from the Main Navigation Menu.

A “Save as PDF” dialog box will be displayed. This allows navigation through the computer to locate the desired patient studies folder. The user can create a destination folder anywhere on the computer or a shared drive. In the example below the folder “Restech Saved PDF” was created in the user’s Documents folder.

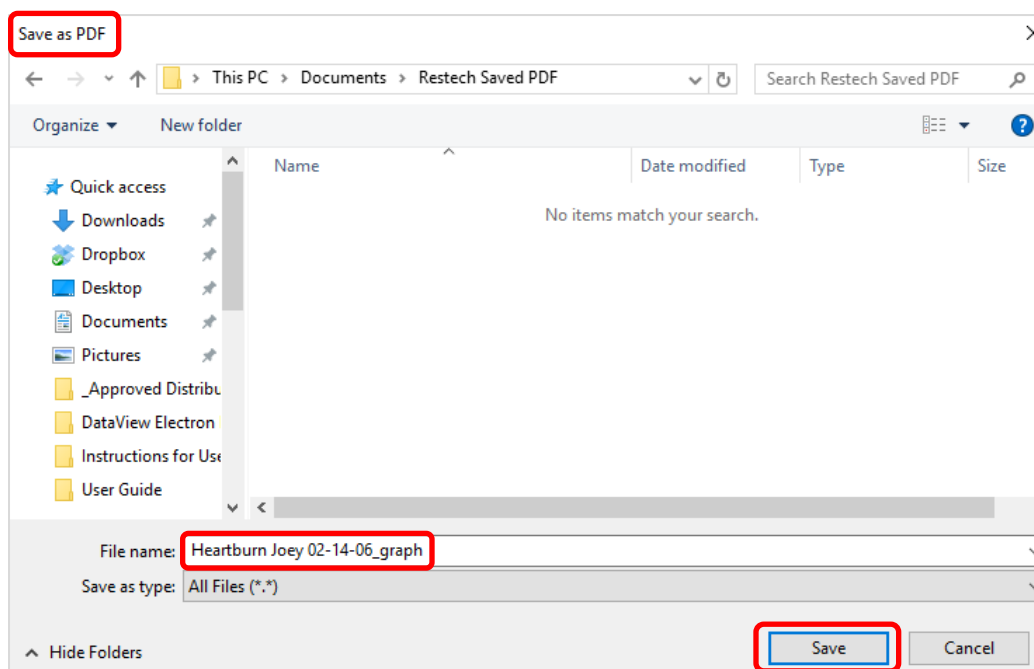


DataView will automatically name the file with last name, first name, study date, and graph.

Select “Save” and the Graph file will be saved in the selected destination folder.


After saving the file as a PDF, the Graph is displayed as a PDF document, using the available PDF program on the computer. DataView automatically selects Landscape Layout for displaying and printing the Graph.

## DataView v4 Software User Guide

The menu items and controls displayed will vary based on the PDF program on the computer.

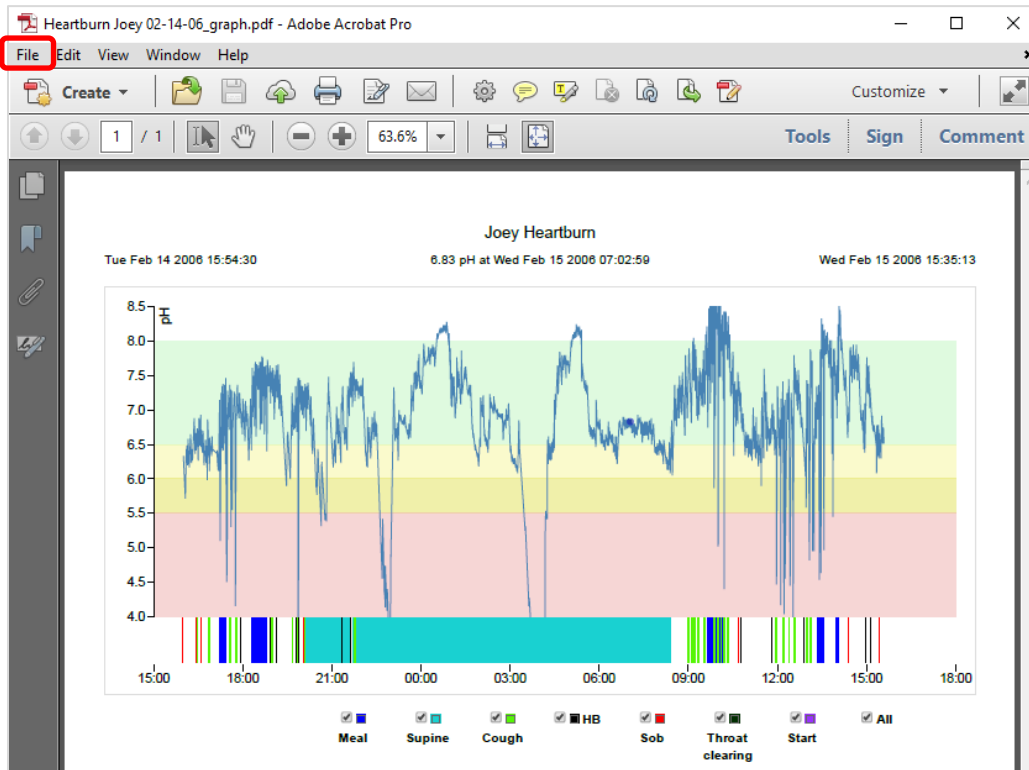


In the opened Graph PDF, select "File" from the menu. Then select "Print" for the desired printer. Follow the Print menu in the PDF program to Print the Graph.

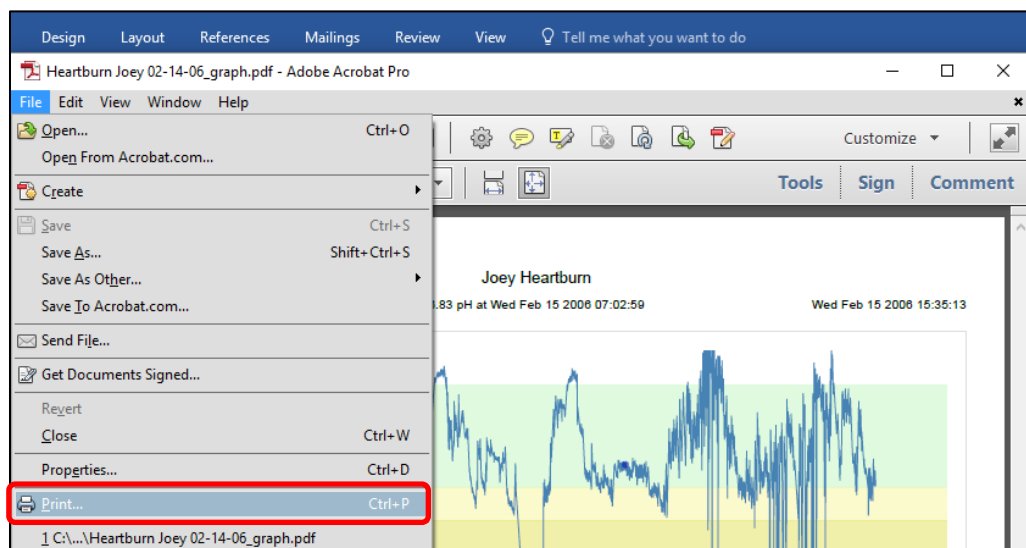

# DataView v4 Software User Guide


## 7.2 Print Report

The "Print" Main Navigation Menu Selection allows the user to save a PDF document of Graphs and Reports.

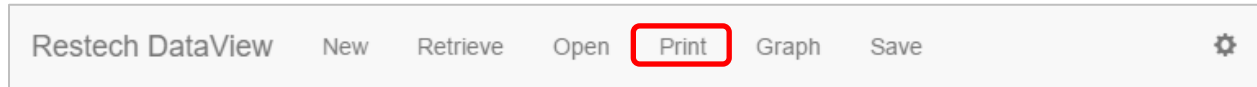

To print a Report from the Report Workspace, select "Print" from the Main Navigation Menu.

A "Save as PDF" dialog box will be displayed. This allows navigation through the computer to locate the desired patient studies folder. The user can create a destination folder anywhere on the computer or a shared drive. In the example below the folder "Restech Saved PDF" was created in the user's Documents folder. The Graph was also saved as a PDF in this folder.

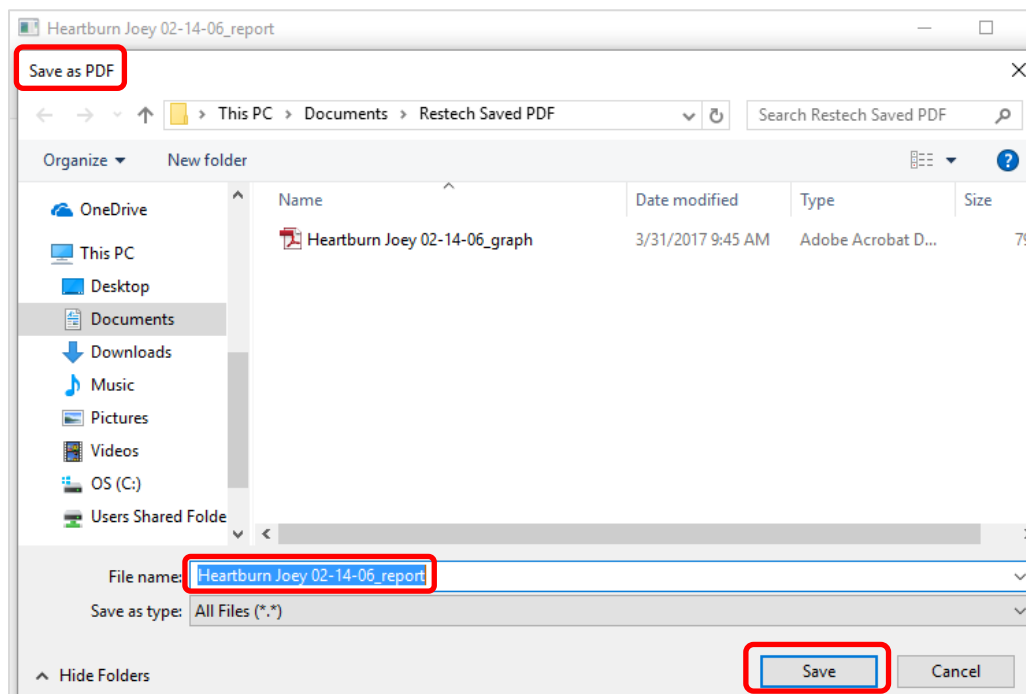


DataView will automatically name the file with last name, first name, study date, and report.

Select "Save" and the Report file will be saved in the selected destination folder.


After saving the file as a PDF, the Report is displayed as a PDF document, using the available PDF program on the computer. DataView automatically selects Portrait Layout for displaying and printing the Report.

# DataView v4 Software User Guide

The menu items and controls displayed will vary based on the PDF program on the computer.

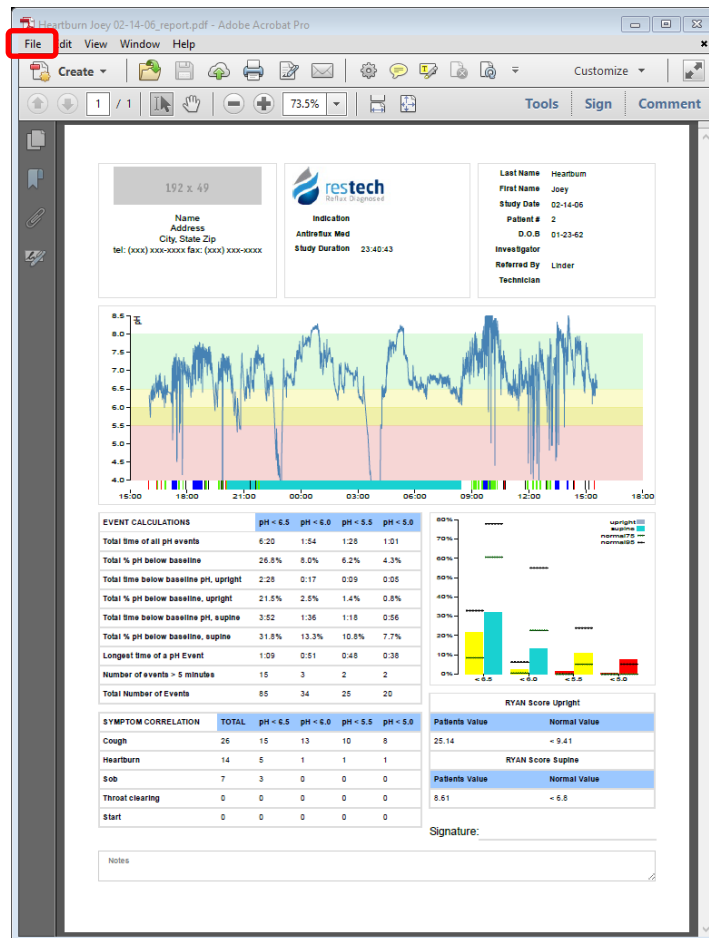

In the opened Report PDF, select “File” from the menu. Then select “Print” for the desired printer. Follow the Print menu in the PDF program to Print the Report.

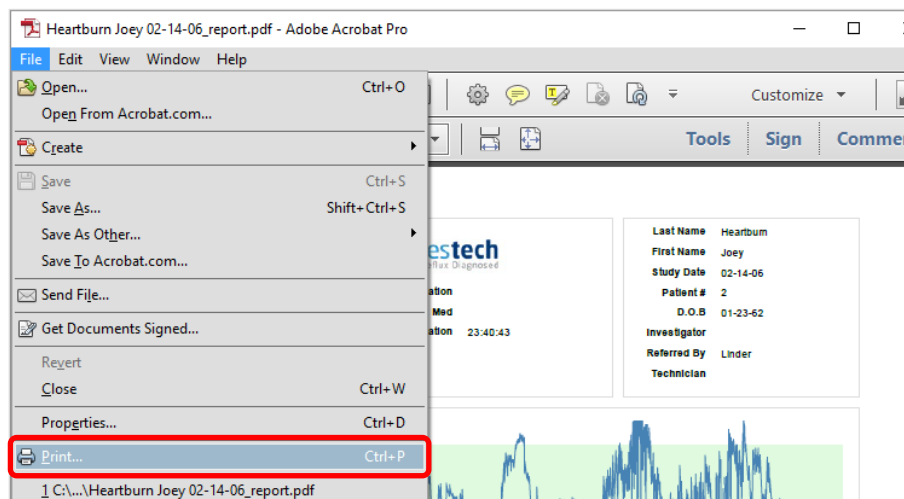

## 8. Main Navigation Menu Selection "Graph"

### 8.1 View and Edit pH Study Graph

The default representation of the patient study files is given in the Graph Workspace. The pH Study Graph displays the pH values versus time. The Event Graph below the pH Graph displays patient symptoms and pH study events.

The "Graph" and "Report" Main Navigation Menu Selections toggle between the Graph and Report Workspaces.




The Patient Name is displayed at the top of the Graph Workspace above the Graph toolbar. The study start time and end time are displayed above the pH Graph.




## DataView v4 Software User Guide

The vertical axis of the pH Study Graph displays pH at 0.5 pH increments and the horizontal axis displays time at 3 hour increments.

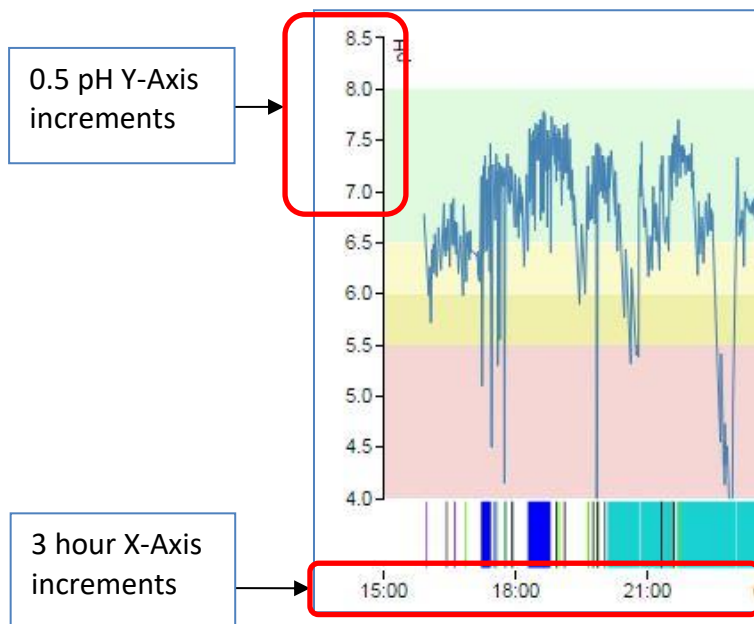

Four horizontal color bands appear on the pH Study Graph for each pH baseline. The colors of these bands may be changed in the "Settings" Main Navigation Menu Selection.

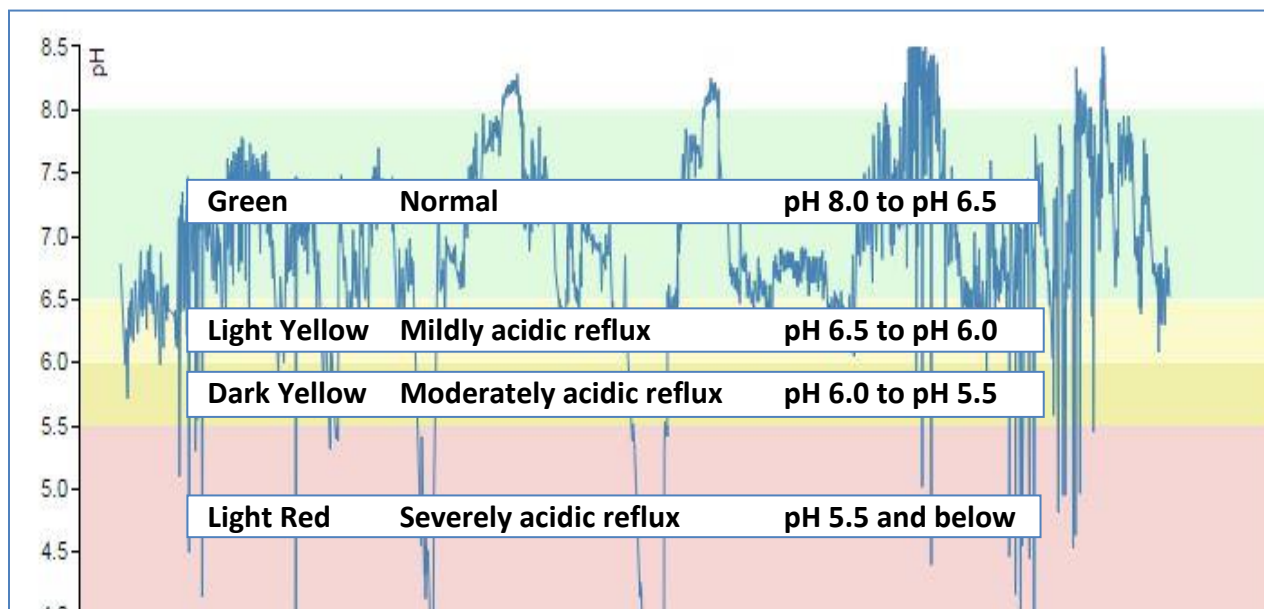


## DataView v4 Software User Guide

### Trackball

The trackball function allows the user to locate specific pH data points or time within the study using the cursor.

Hover the mouse cursor over the pH data of the pH Graph and move horizontally along the points of the study. The trackball follows the movement of the mouse cursor along the pH Graph.

The trackball pH and time are displayed above the pH Graph in the area below the Graph toolbar. The pH value and time on the trackball display change to match the pH, time, and date of each point on the pH Graph.

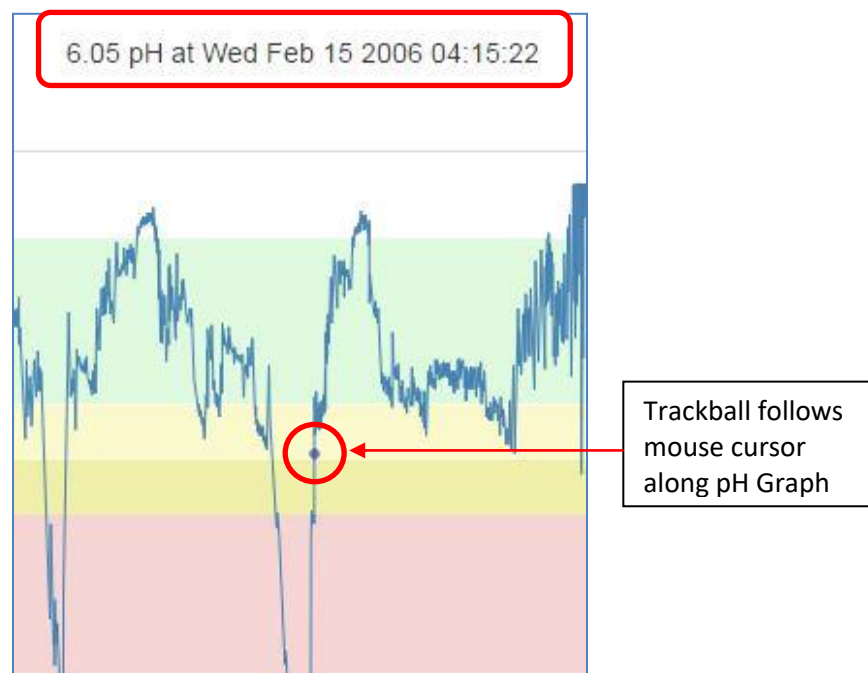

### Missing time periods and dry out events

If the data transmission is interrupted during the study, the pH Graph will display a straight line connecting the last transmitted pH value and the next transmitted pH value.

If the pH Probe dries out during the study, the pH Graph will display a straight line connecting the last transmitted pH value and the next transmitted pH value.


If the pH Probe is removed from the patient and dries out before the pH study ends, the pH Graph will display a straight line from the last transmitted pH value to the end of the study.




## DataView v4 Software User Guide

### 8.2 View and Edit Event Graph


The Event Graph displays events marked by patient button presses on the Dx-Recorder during the pH study. You may visually correlate patient symptoms, meals, and supine periods with changes in pH. The pH Graph and Event Graph share the X-axis (time).

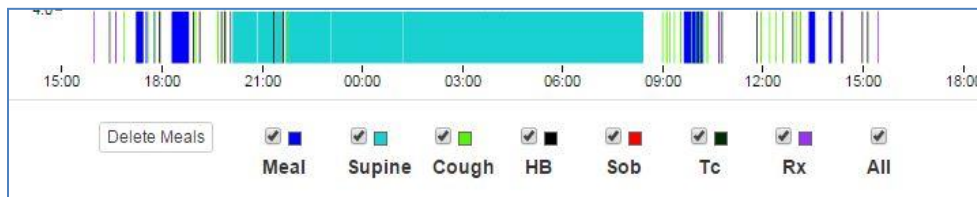


Each checkbox allows the user to show or hide the corresponding symptom. Unchecking the meals checkbox removes the blue meal markers from the Event Graph.

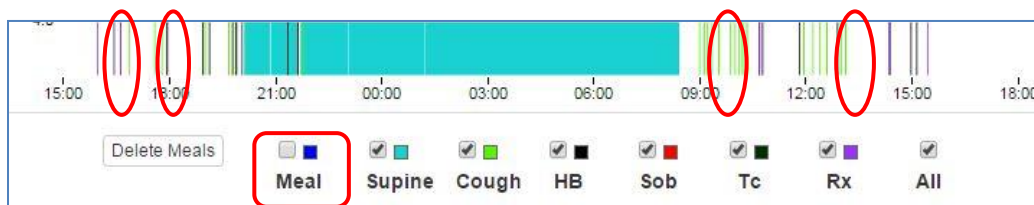


Unchecking the supine checkbox removes the teal supine markers from the Event Graph.

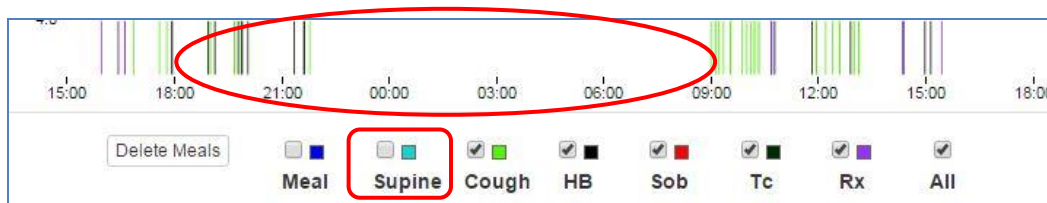

Unchecking the symptom checkboxes removes the button press symptom markers from the Event Graph.

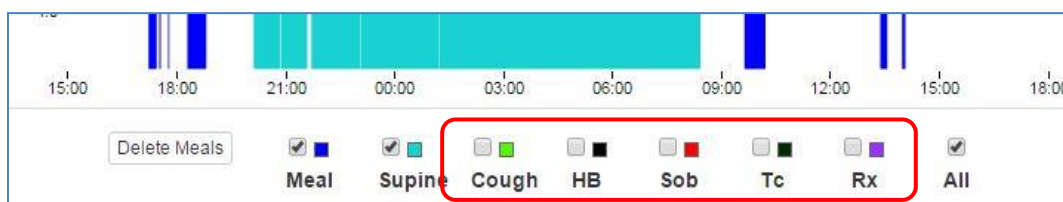


Unchecking the All checkbox removes all event and symptom markers from the Event Graph.

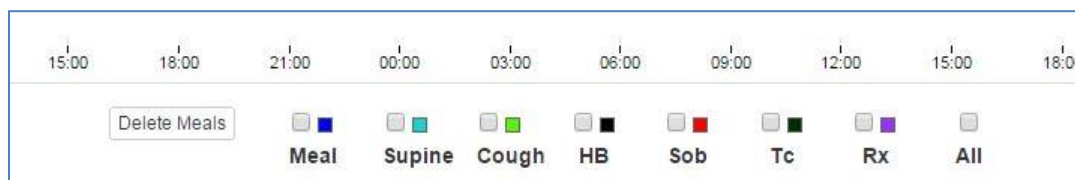

## DataView v4 Software User Guide

### 8.3 Graph Toolbar


The DataView Software contains tools located in the Graph toolbar to precisely edit the pH Graph to match the written patient diary entries.

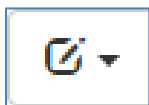


Add events and symptoms in a selected range



Remove events and symptoms in a selected range — does not remove pH data

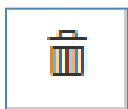

Delete pH data, events, and symptoms in a selected range — removes pH data

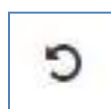


Undo the last undoable command: erases the last change done to the Graph reverting it to the previous state



Zoom in on pH data in selected range



Zoom out to original size of Graph

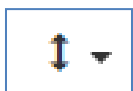

Adjust the Y-axis pH scale


**NOTE: CHECK the patient's diary entries before editing the pH Graph.  
COMPARE THE EVENT GRAPH TO THE WRITTEN PATIENT DIARY ENTRIES.**

## DataView v4 Software User Guide


### 8.3.1 Add Event and Symptom Information

The Add button in the Graph toolbar allows the user to add Events and Symptoms to the Graph.


Hover the mouse over the Add button and the tooltip displays "Add period to pH data in a range."

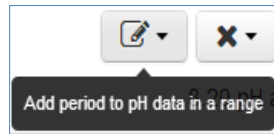

Select a time range on the pH Graph using the click and drag function. Click the Add button. The Add button displays a dropdown menu containing a list of available Events and Symptoms. The last three Custom Symptoms will vary depending on the symptoms input for each patient study.

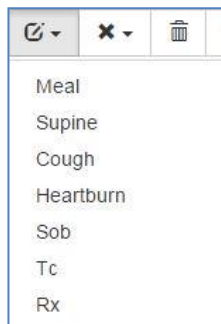

If a meal is added to a selected range, a new meal marker will appear in the Event Graph.

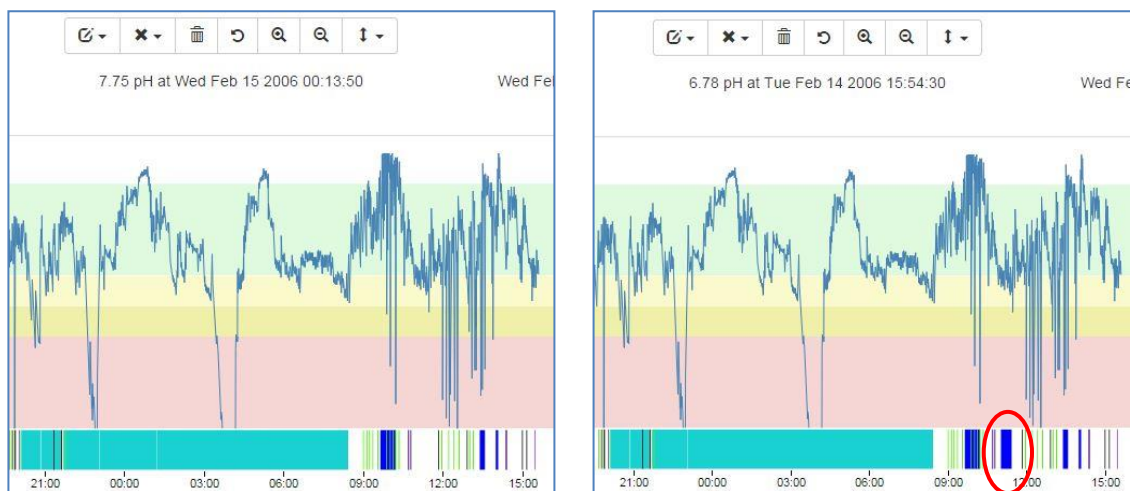


## DataView v4 Software User Guide


### 8.3.2 Remove Event and Symptom Information

The Remove button in the Graph toolbar allows the user to remove Events and Symptoms on the Graph. The data is not removed from event calculations.


Hover the mouse over the Remove button and the tooltip displays "Clear periods from pH data in a range."

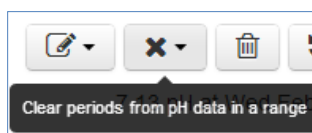

Select a time range on the pH Graph using the click and drag function. Click the Remove button. The Remove button displays a dropdown menu containing a list of available Events and Symptoms. The last three Custom Symptoms will vary depending on the symptoms input for each patient study. Clicking "All" removes all Events and Symptoms.

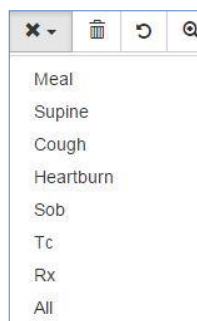

If a meal is removed in a selected range, a meal marker will disappear in the Event Graph.

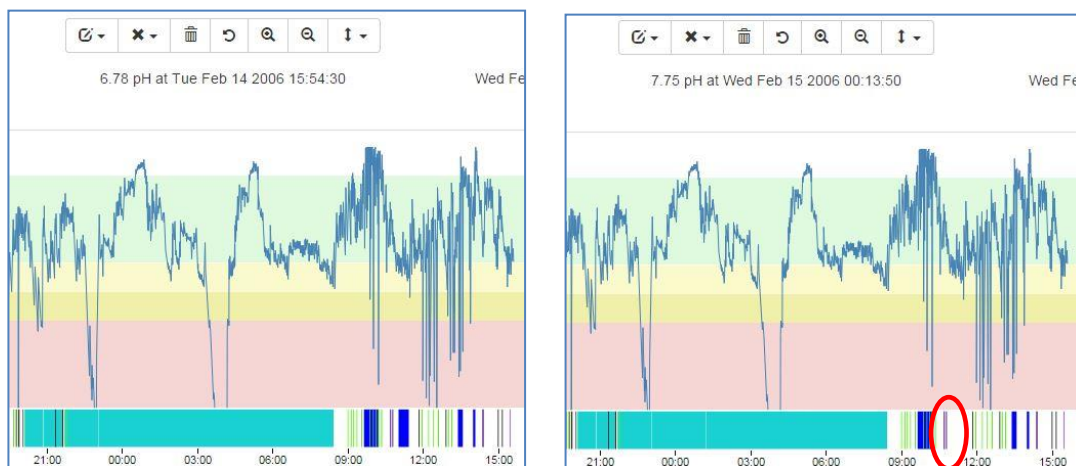


## DataView v4 Software User Guide

### 8.3.3 Delete Event Information


The Delete button in the Graph toolbar allows the user to delete pH data, Events, and Symptoms on the Graph. The data is removed from event calculations.

Hover the mouse over the Delete button and the tooltip displays "Delete pH data and periods in a range."

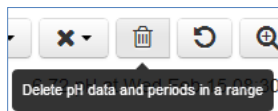


Select a time range on the pH Graph using the click and drag function. Click the Delete button. The Delete button completely removes pH data, Events, and Symptoms in the selected range.

If a time range that includes a meal period is deleted, all of the pH data, including the meal, is completely deleted.

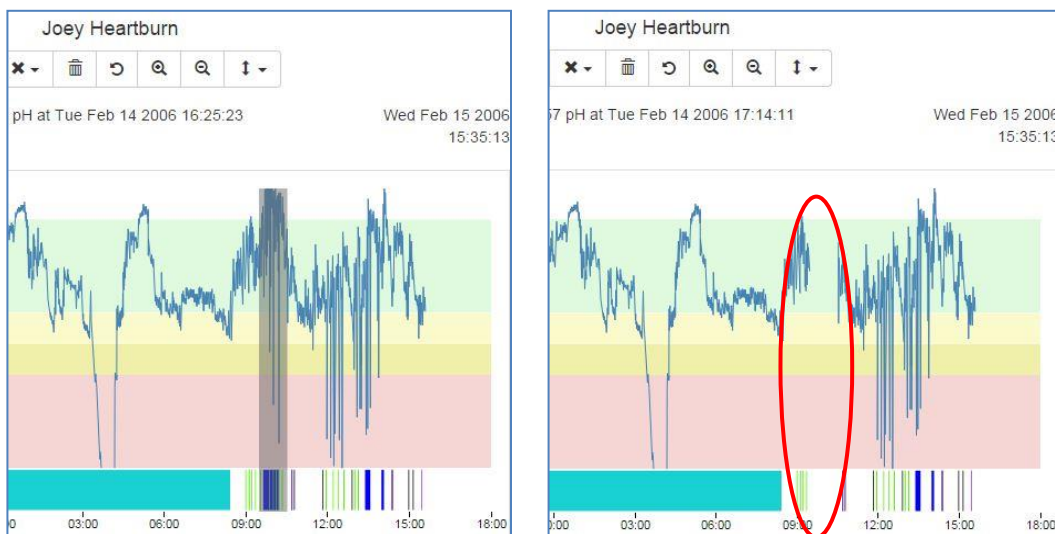

### 8.3.4 Undo Changes

The Undo button in the Graph toolbar allows the user to undo commands in the Graph function. Hover the mouse over the Undo button and the tooltip displays "Undo the last undoable command."

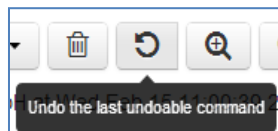


Click the Undo button. The Undo button reverts to the last change made. Click the Undo button multiple times in order to undo all the changes made during the current session.


## DataView v4 Software User Guide

### 8.3.5 Zoom In

The Zoom In button in the Graph toolbar allows the user to view a magnified section of the Graph by expanding the X-axis. Hover the mouse over the Zoom In button and the tooltip displays "Zoom in on pH data in a selected range."

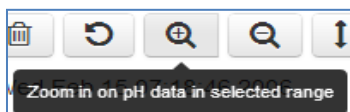

Select a time range on the pH Graph using the click and drag function. Click the Zoom In button. The Zoom In button magnifies the pH Graph data of the selected range.

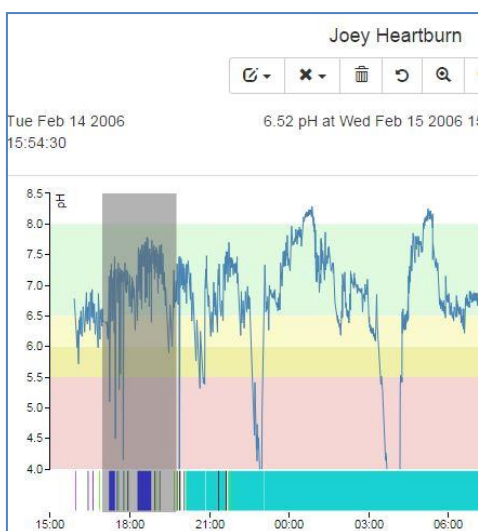


Select a time range



Zoomed in view of selected range

### 8.3.6 Zoom Out

The Zoom Out button in the Graph toolbar allows the user to restore the graph display to the original time interval of the Graph. Hover the mouse over the Zoom Out button and the tooltip displays "Zoom out."

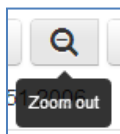

Click the Zoom Out button. The Zoom Out button returns the pH Graph to its original size.

## DataView v4 Software User Guide




Zoomed in view




Zoom out to return to original size

### 8.3.7 Y-Axis (pH Scale) Adjustment

The Y-Axis Adjustment button in the Graph toolbar allows the user to adjust the Y-Axis pH scale of the Graph. Hover the mouse over the Y-Axis Adjustment button and the tooltip displays "Adjust pH scale."

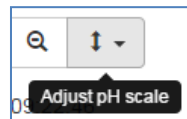

Click the Y-Axis Adjustment button. The Y-Axis Adjustment button displays a dropdown menu containing a list of available pH ranges: 4 to 8.5, 0 to 10, and Custom.

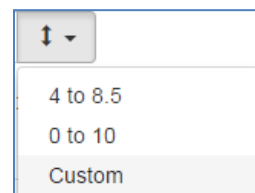


Select the 0 to 10 pH setting from the dropdown menu. The Graph Y-Axis changes to 0 to 10 pH.






## DataView v4 Software User Guide


Select the Custom pH setting from the dropdown menu. A Y-Axis Adjustment dialog box is displayed for entering the custom Top pH value. The number 9 is entered as an example here. Click "OK" to continue.

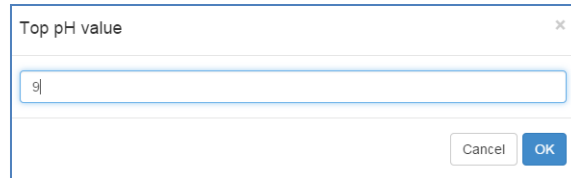A dialog box titled "Top pH value" with a close button (X) in the top right corner. It contains a text input field with the number "9" entered. At the bottom right, there are two buttons: "Cancel" and "OK".


A Y-Axis Adjustment dialog box is displayed for entering the custom Bottom pH value. The number 1 is entered as an example here. Click "OK" to continue.

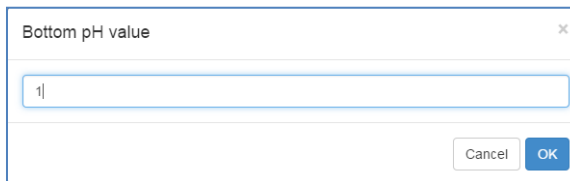A dialog box titled "Bottom pH value" with a close button (X) in the top right corner. It contains a text input field with the number "1" entered. At the bottom right, there are two buttons: "Cancel" and "OK".

The Y-Axis Adjustment dialog box closes and the pH Graph changes to the new Y-Axis scale.

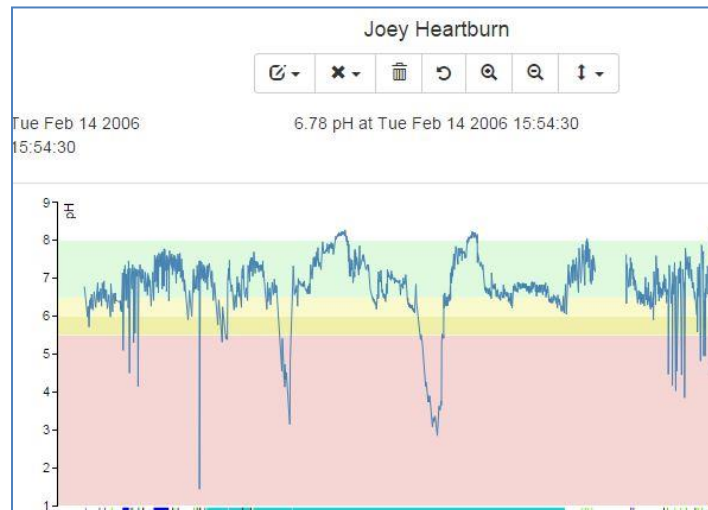

To return to the default Y-Axis, click the Y-Axis Adjustment button and select 4 to 8.5 in the dropdown menu.

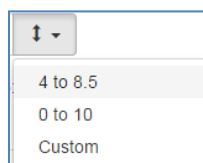A dropdown menu with a small icon of a vertical axis with an arrow pointing up and down. The menu is open, showing three options: "4 to 8.5", "0 to 10", and "Custom".




## DataView v4 Software User Guide


### 8.3.8 Delete Meals

The Delete Meals button in the Graph Workspace allows the user to delete pH data associated with all meals on the pH Graph.


Clicking the Delete Meals button changes the calculations for the Event Calculations Table, but *does not remove the meal markers* from the Event Graph and *does not remove pH data* from the pH Graph.

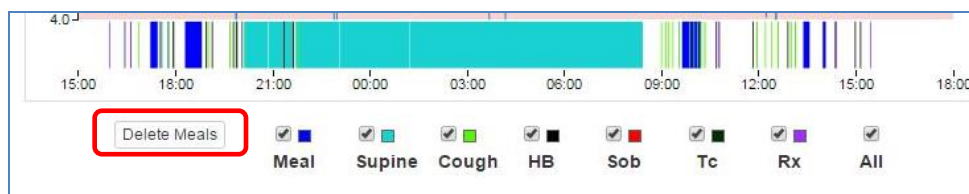

Click the Delete Meals button. The 'Delete how many minutes before/after each meal period?' dialog box is displayed. The default number is 5 minutes, but the number may be changed in the selection field. Click "OK" to continue.

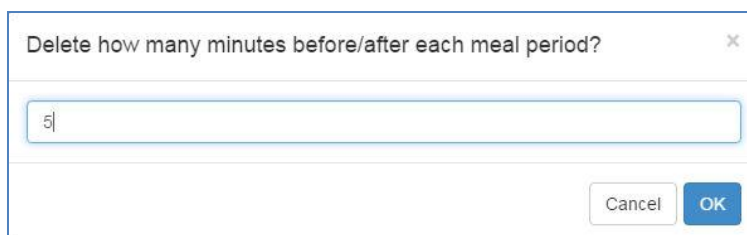

The Delete Meals button now displays "Undelete Meals" and the Event Graph still shows the meal events.

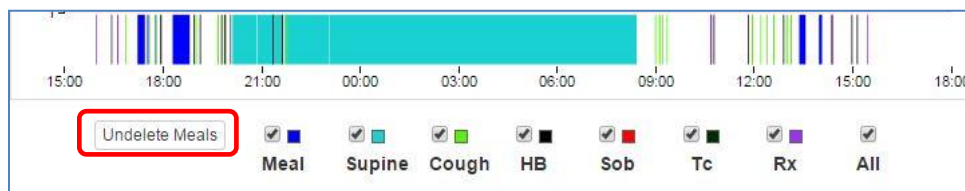


To confirm that meals have been deleted from the calculations, click Report on the Main Navigation Menu and locate the center header. The Restech Information Block displays 'Meals Removed from Calculations.'




## 9. Main Navigation Menu Selection “Report”

### 9.1 View and Edit Reports

The "Report" Main Navigation Menu Selection displays a representation of the patient study file in the Report Workspace.

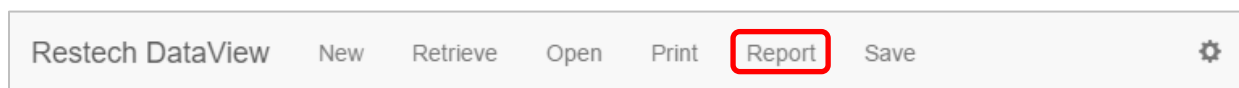

The Report Workspace displays the Site Information, Restech Information, Patient Study Information, pH Study Graph, Acid Exposure Graph, Event Calculation Table, Symptom Calculation Table, RYAN Score, Notes field, and Signature field.

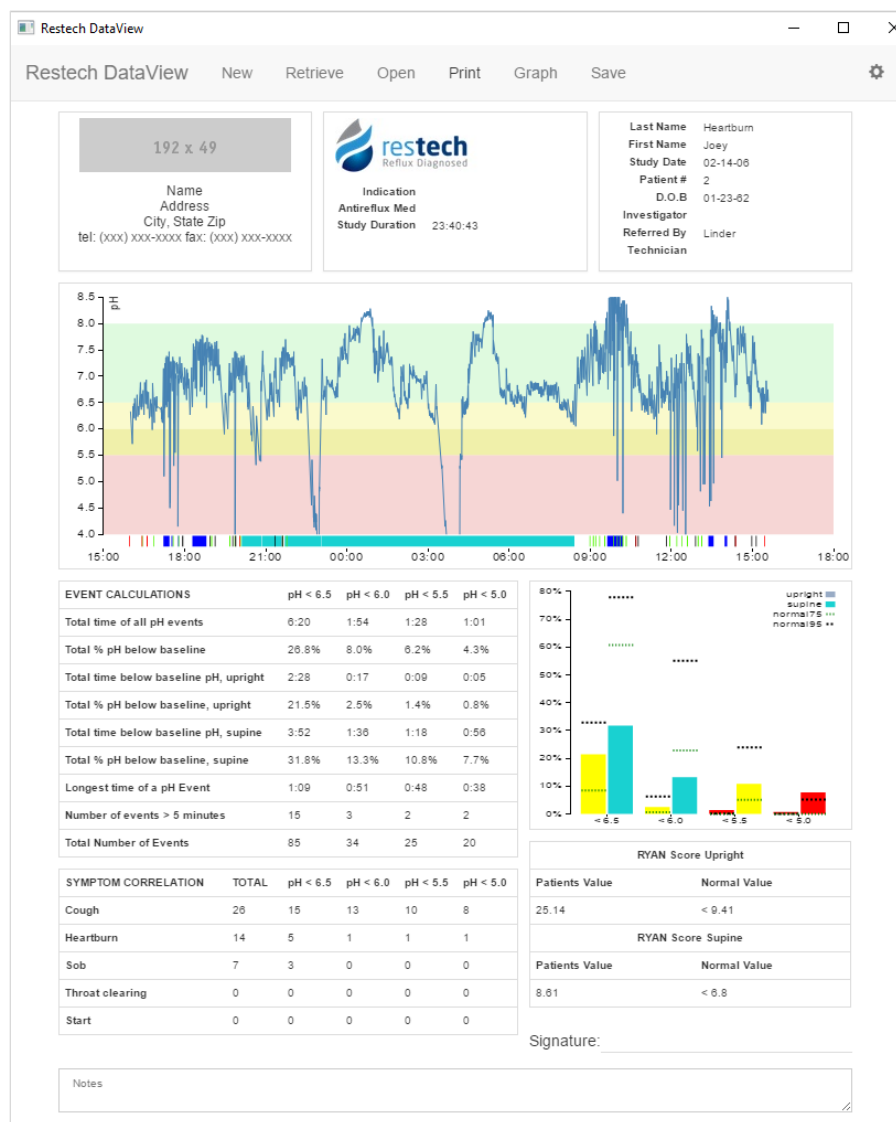

## DataView v4 Software User Guide

### 9.1.1 Site Information Block

The Site Information Block displays the site logo image and header in the upper left corner of the Report Workspace.

The user has the option to upload an image and edit site information in Settings. See Section 11.1.1 for instructions to edit and save site information.

Once an image is uploaded and edited using site information in Settings, the changes will be saved on DataView until the software is uninstalled. The user may edit and overwrite site information saved in Settings.



### 9.1.2 Restech Information Block

The Restech Information Block displays the Restech logo image and header, patient indications and anti-reflux medication, and study duration in the upper middle section of the Report Workspace.



If meals have been deleted from the calculations, the Restech Information Block displays 'Meals Removed from Calculations.'



### 9.1.3 Patient Study Information Block

The Patient Study Information Block displays the Patient Details and Study Details inputs in the upper right corner of the Report Workspace.

Information in the Patient Study Information Block is saved from the "New" Main Navigation Menu Selection using the Patient Details Tab and Study Details Tab.


Last Name	Heartburn
First Name	Joey
Study Date	02-14-08
Patient #	2
D.O.B	01-23-82
Investigator	
Referred By	Linder
Technician	

### 9.1.4 pH Study Graph

The pH Study Graph displays a reduced scale replica of the pH Study Graph and Event Graph displayed on the Graph Workspace.


Unlike the pH Study Graph and Event Graph on the Graph Workspace, the pH Study Graph on the Report Workspace is not editable. Changes made on the Graph Workspace will automatically appear on the Report Workspace.

The Y-Axis of the pH Study Graph on the Graph Workspace determines the Y-Axis of the pH Study Graph on the Report Workspace.

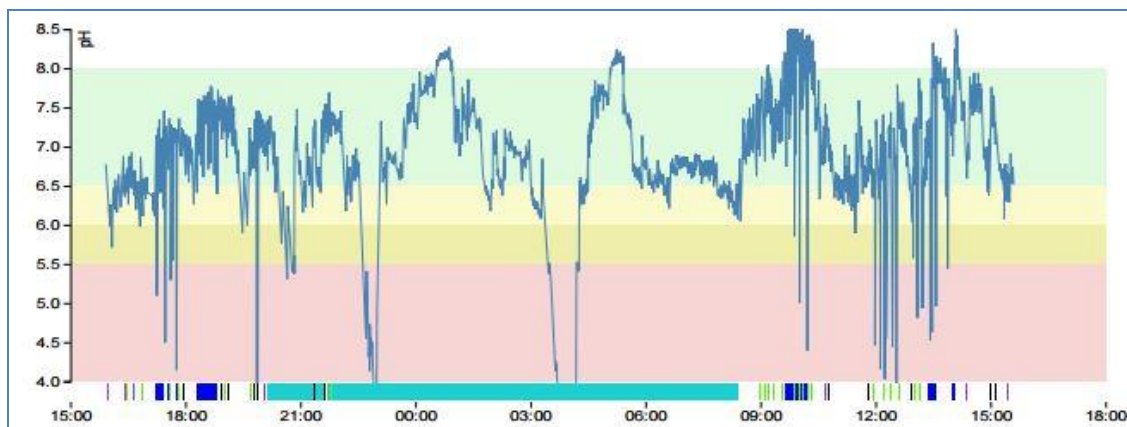

## 9.1.5 Acid Exposure Graph

The Acid Exposure Graph displays the “Total % below baseline, upright” and “Total % below baseline, supine” in bar graph format. These values are in the Event Calculations Table (see Section 9.1.6). Upright (left bar) and supine (right bar) are displayed at the four baseline pH ranges presented in the Event Calculations Table.

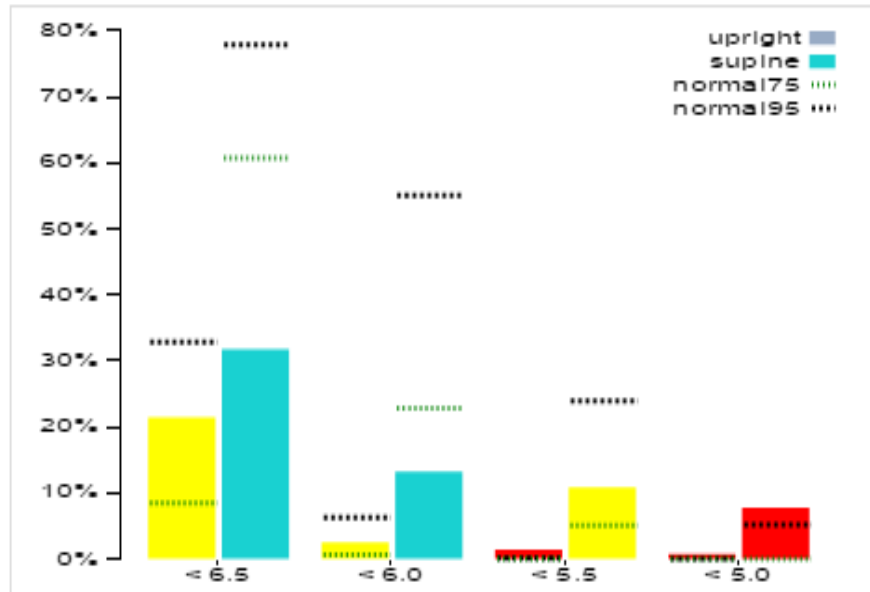

Example Acid Exposure Graph

The Y-Axis of the Acid Exposure Graph is fixed at 80% for all patients.  
The X-Axis of the Acid Exposure Graph is fixed for each pH baseline.

### Normal Value Markers

75<sup>th</sup> percentile = Thin horizontal dashed lines are marked for each pH baseline

95<sup>th</sup> percentile = Bold horizontal dashed lines are marked for each pH baseline

### Color Key Upright

Gray = normal; patient's total % below baseline does not exceed the acid normal values.

Yellow > 75<sup>th</sup> percentile; patient's total % below baseline exceeds 75<sup>th</sup> percentile normal value.

Red > 95<sup>th</sup> percentile; patient's total % below baseline exceeds 95<sup>th</sup> percentile normal value.

### Color Key Supine

Teal = normal; patient's total % below baseline does not exceed the acid normal values.

Yellow > 75<sup>th</sup> percentile; patient's total % below baseline exceeds 75<sup>th</sup> percentile normal value.

Red > 95<sup>th</sup> percentile; patient's total % below baseline exceeds 95<sup>th</sup> percentile normal value.

More information about the percentile ranges can be found in the Normative DataChart in the Settings → Help menu.

## DataView v4 Software User Guide

### 9.1.6 Event Calculations Table

The Event Calculations Table automatically calculates and displays events on the Report Workspace for each baseline in both the upright and supine position.

EVENT CALCULATIONS	pH < 6.5	pH < 6.0	pH < 5.5	pH < 5.0
Total time of all pH events (hh:mm)	6:20	1:54	1:28	1:01
Total % pH below baseline pH	26.8%	8.0%	6.2%	4.3%
Total time below baseline pH, upright	2:28	0:17	0:09	0:05
Total % pH below baseline, upright	21.5%	2.5%	1.4%	0.8%
Total time below baseline pH, supine	3:52	1:36	1:18	0:56
Total % below baseline pH, supine	31.8%	13.3%	10.8%	7.7%
Longest time of a pH event	1:09	0:51	0:48	0:38
Number of events > 5 minutes	15	3	2	2
Total Number of Events	85	34	25	20

Event calculations for each baseline are not cumulative. They are standalone calculations at each pH range for percent, time, and number of events.

Example: If the patient's pH goes below 6.0 pH and then continues further below 5.5 pH and then below 5.0 pH, that pH "event" will be counted as time below 6.0 pH, then as time below 5.5 pH, and then as time below 5.0 pH for the period of time the pH was below each of those levels. When the pH increases to above a pH level, it is no longer counted for that pH level.

Total time and total % pH below baseline in the upright position are calculated from pH data in the upright position only. The percent of time below baseline upright is divided by the amount of upright time in the pH study.

Similarly, total time and total % pH below baseline in the supine position are calculated from pH data in the supine position only. The percent of time below baseline supine is divided by the amount of supine time in the pH study.

Clicking the 'Delete Meals' button on the Graph Workspace deletes the pH data during the meal including the selected pre-and post-prandial range, changing the affected calculations presented in the Event Calculation Table.

## DataView v4 Software User Guide

### 9.1.7 Symptom Correlation Table

The Symptom Correlation Table automatically calculates and displays symptoms on the Report Workspace for each baseline in both the upright and supine position.

SYMPTOM CORRELATION	TOTAL	pH < 6.5	pH < 6.0	pH < 5.5	pH < 5.0
cough	26	15	13	10	8
heartburn	14	5	1	1	1
Sob	0	0	0	0	0
Throat clearing	0	0	0	0	0
Start	7	3	0	0	0

The “Total” column displays total button presses for each symptom. Values displayed in the table for each pH level indicate how many times the patient pressed that symptom button while in that pH range. In the example above, there were 13 cough button presses during the period when the patient’s pH was below 6.0 pH. There were 10 cough button presses when the patient’s pH was below 5.5 pH, and 8 cough button presses when the patient’s pH was below 5.0 pH. There were 26 cough button presses during the pH study.

Clicking the 'Delete Meals' button on the Graph Workspace deletes the pH data during the meal including the selected pre-and post-prandial range, changing the affected calculations presented in the Symptom Correlation Table.

### 9.1.8 RYAN Score Table

The RYAN Score Table displays a score for both the upright and supine positions. The RYAN Score published values are displayed as normal values.



RYAN Score Upright	
Patient's Value	Normal Value
25.14	< 9.41
RYAN Score Supine	
Patient's Value	Normal Value
8.61	< 6.8

More information about the RYAN Score can be found in the DeMeester RYAN Score article in the Settings → Help menu.

## DataView v4 Software User Guide

### 9.1.9 Notes


The Notes field displays an editable text box that allows the user to directly type in notes. The notes will be saved and displayed on the printed Report.

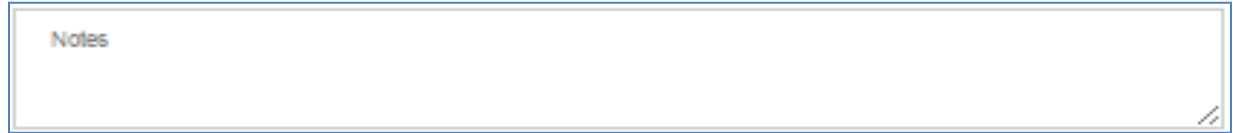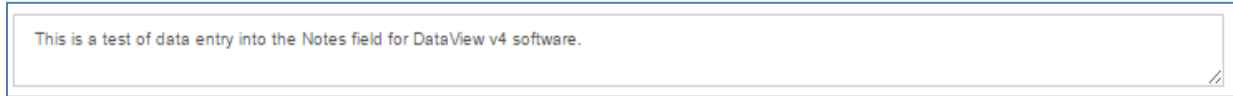A screenshot of a rectangular text box with a thin blue border. The word "Notes" is written in a small, light blue font in the top-left corner. In the bottom-right corner, there is a small icon consisting of two parallel diagonal lines.A screenshot of a rectangular text box with a thin blue border, similar to the one above. It contains the text "This is a test of data entry into the Notes field for DataView v4 software." in a small, light blue font. A small icon of two parallel diagonal lines is visible in the bottom-right corner.

### 9.1.10 Signature

The Signature field allows the user to manually enter a signature on the printed Report. The signature field will appear blank on the Report Workspace.

A signature cannot be directly input from the Report Workspace.

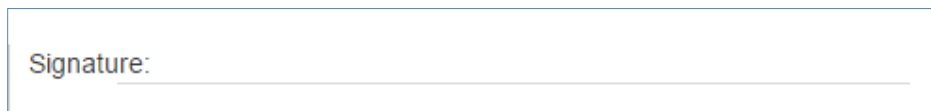
A screenshot of a rectangular text box with a thin blue border. The word "Signature:" is written in a small, light blue font on the left side. A horizontal line is drawn across the box, starting from the text and extending to the right edge.




### 10. Main Navigation Menu Selection "Save"

The "Save" Main Navigation Menu Selection allows the user to save a modified version of the patient study after editing the Graph.


Note: A patient study retrieved from a RESTECH formatted SD card will automatically be saved. Do not use "Save" when retrieving a patient study from the RESTECH SD card.

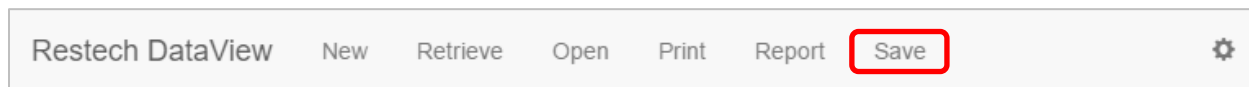

Hover the mouse over the "Save" Main Navigation Menu Selection and the tooltip displays "Save study to disk."

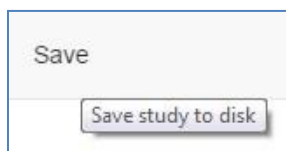

If modifications have been made to a pH study, click "Save" in the Main Navigation Menu. The 'Please wait Saving' box is displayed.

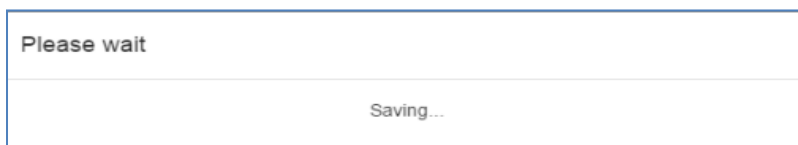

The modified study will be saved using the patient's last name, first name, study date, followed by 'modified'.

The original unmodified patient study will remain unchanged in Restech DataView data folder.


Once the modified study is saved to the computer, the 'Please wait Saving' box disappears and the DataView Workspace is displayed.

To verify that a modified study has been saved, click "Open" in the Main Navigation Menu and the modified patient study will be displayed in the Open Study dialog box.


## 11. Main Navigation Menu Selection “Settings”

### 11.1 Settings

The "Settings" Main Navigation Menu Selection provides users with submenu items.

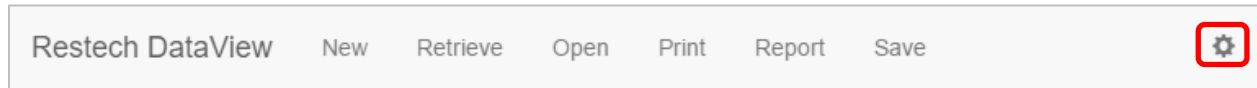

Click "Settings" in the Main Navigation Menu. The settings menu displays a dropdown menu with the submenu items: 'Settings' and 'Help.'

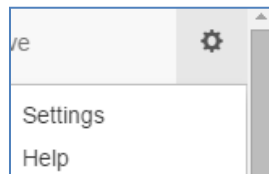

#### 11.1.1 General Tab

Click on the 'Settings' submenu item to open the Settings dialog box. The first tab displays General input fields.

Click in the entry fields to input Site Information data according to the available input fields.

The first field allows the user to upload an image logo.

## DataView v4 Software User Guide

Information inserted into the General Tab will appear on the Site Information Block in the upper left corner of the Report Workspace.



A rectangular placeholder box with a blue border. Inside, there is a grey rectangle at the top labeled "192 x 49". Below it, the text "Name", "Address", "City, State Zip", "tel: (xxx) xxx-xxxx fax: (xxx) xxx-xxxx" is displayed in a centered, sans-serif font.


Once a user has uploaded an image and edited Site Information in Settings, the changes will be saved in DataView until the software is uninstalled. The user may edit and overwrite Site Information by returning to the General Tab.

### 11.1.2 Normal Acid Exposure

Click on the 'Settings' submenu item to open the Settings dialog box. The second tab displays Acid Normals. These are the pharyngeal pH values for normal patients at the 95<sup>th</sup> and 75<sup>th</sup> percentile in upright and supine.

The values are extracted from the publication “A New Technique for Measurement of Pharyngeal pH: Normal Values and Discriminating pH Threshold” by Ayazi, et al, available in the DataView Help menu link “DeMeester RYAN Score article.”

The columns in the Acid Normals tab represent pH values < 6.5, < 6.0, < 5.5, and < 5.0 as displayed in the Event Calculations Table. These values cannot be edited by the user.

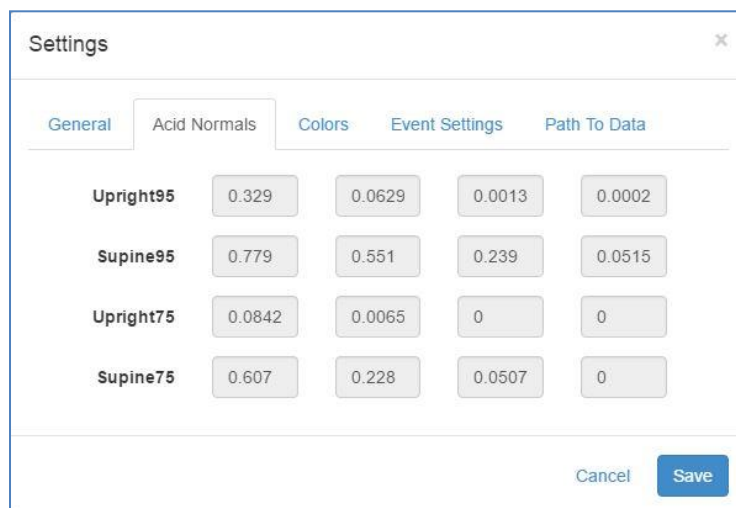

The Settings dialog box is shown with the 'Acid Normals' tab selected. It contains a table of pH values for four scenarios: Upright95, Supine95, Upright75, and Supine75. Each scenario has four corresponding input fields with numerical values. At the bottom right, there are 'Cancel' and 'Save' buttons.

	< 6.5	< 6.0	< 5.5	< 5.0
Upright95	0.329	0.0629	0.0013	0.0002
Supine95	0.779	0.551	0.239	0.0515
Upright75	0.0842	0.0065	0	0
Supine75	0.607	0.228	0.0507	0


## DataView v4 Software User Guide

### 11.1.3 Colors Tab

Click on the 'Settings' submenu item to open the Settings dialog box. The third tab displays “Colors” input fields.

These values may be edited by the user by clicking in the desired field and typing in new values.

Once the Colors have been edited, the changes will be saved on DataView until the software is uninstalled. The user may edit and overwrite the Colors by returning to the Colors Tab. Return to the default colors by selecting Load Default Colors.

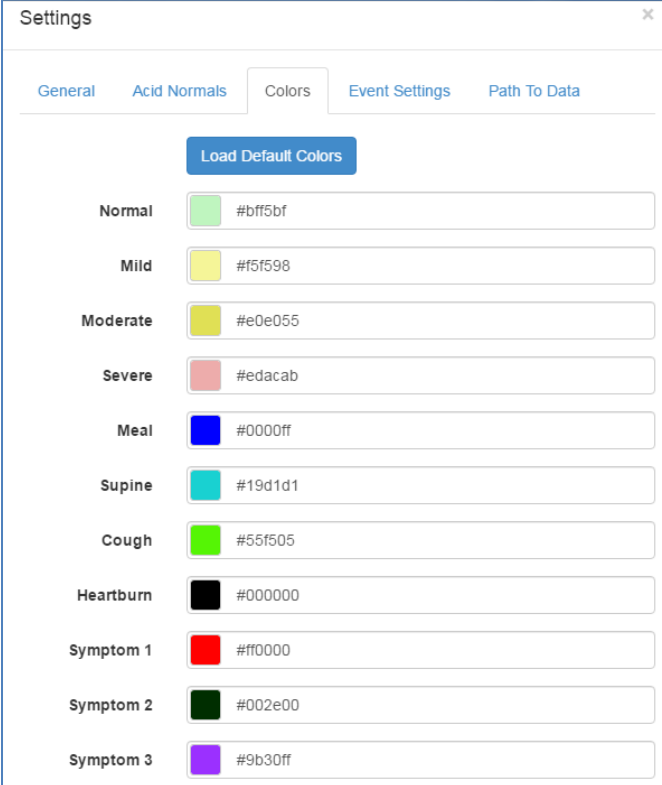

The screenshot shows the 'Settings' dialog box with the 'Colors' tab selected. The dialog has a title bar with a close button. Below the title bar are five tabs: 'General', 'Acid Normals', 'Colors', 'Event Settings', and 'Path To Data'. The 'Colors' tab is active, showing a 'Load Default Colors' button at the top. Below this button are ten rows, each with a label, a color swatch, and a text input field containing a hex color code. The labels are: Normal, Mild, Moderate, Severe, Meal, Supine, Cough, Heartburn, Symptom 1, Symptom 2, and Symptom 3. The color swatches are: light green, yellow, olive green, pink, blue, cyan, bright green, black, red, dark green, and purple. The hex codes are: #bfff5b, #f5f598, #e0e055, #edacab, #0000ff, #19d1d1, #55f505, #000000, #ff0000, #002e00, and #9b30ff.


Label	Color Swatch	Hex Code
Normal	Light Green	#bfff5b
Mild	Yellow	#f5f598
Moderate	Olive Green	#e0e055
Severe	Pink	#edacab
Meal	Blue	#0000ff
Supine	Cyan	#19d1d1
Cough	Bright Green	#55f505
Heartburn	Black	#000000
Symptom 1	Red	#ff0000
Symptom 2	Dark Green	#002e00
Symptom 3	Purple	#9b30ff

These Colors represent the pH Graph, event markers, and symptom markers on both the Graph Workspace and the Report Workspace.

## DataView v4 Software User Guide

### 11.1.4 Event Settings


The fourth tab in the 'Settings' submenu displays Event Settings. When a pH event is above a baseline pH value by 0.1 pH an event timer called the minimum event time “Min Event Time” counts for 10 seconds before ending the pH event. The values cannot be changed by the user.

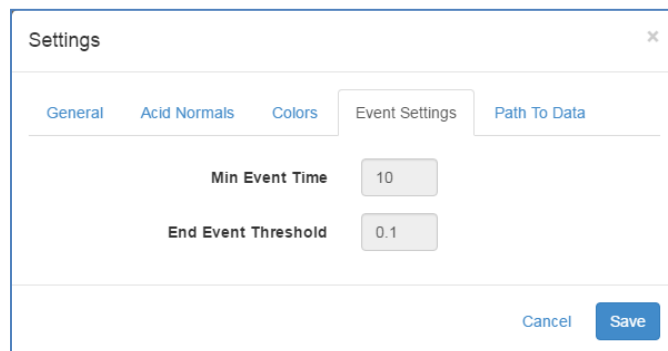

The screenshot shows the 'Settings' dialog box with the 'Event Settings' tab selected. The 'Min Event Time' is set to 10 and the 'End Event Threshold' is set to 0.1. There are 'Cancel' and 'Save' buttons at the bottom right.

### 11.1.5 Path to Data


The fifth tab in the 'Settings' submenu displays Path to Data. This is the location where DataView saves the pH studies retrieved from data cards. The location cannot be edited by the user. Contact Restech Technical Support to modify the Path to Data.

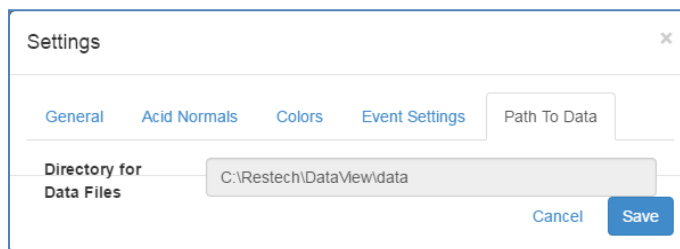

The screenshot shows the 'Settings' dialog box with the 'Path To Data' tab selected. The 'Directory for Data Files' is set to 'C:\Restech\DataView\data'. There are 'Cancel' and 'Save' buttons at the bottom right.

## 11.2 Help

The 'Help' submenu item displays the DataView About, License Agreement, User Guide, DeMeester RYAN Score article, and Normative Data Chart.

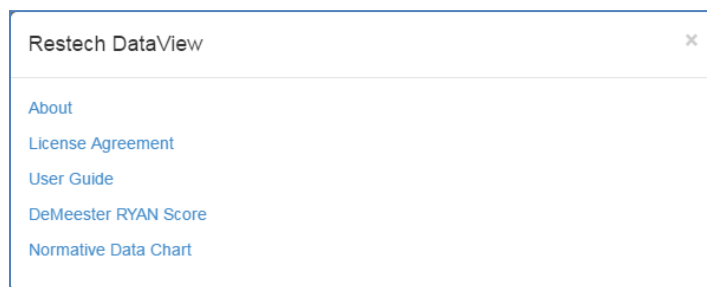

The screenshot shows the 'Restech DataView' help menu with the following options: About, License Agreement, User Guide, DeMeester RYAN Score, and Normative Data Chart.

Click anywhere outside of the menu field to exit the 'Help' submenu.

### 12. Uninstalling DataView v4 Software

The DataView software may be uninstalled to remove DataView from the computer. Uninstall will remove the DataView program files and registry entries.

1	Navigate to Control Panel from the Start Menu and open Programs and Features.
2	Scroll down to find and select the Restech DataView program. Click Uninstall and the Restech DataView software will be uninstalled from the computer.
3	If any DataView shortcuts were pinned to the taskbar or remain after removing DataView, right click on the shortcut icon and click Unpin or Delete to remove it.

#### Trademarks

---

*Restech and DataView are trademarks of Respiratory Technology Corporation.*

*Pentium is a registered trademark of Intel Corporation in the United States and other countries.*

*Windows is a registered trademark of Microsoft Corporation in the United States and other countries.*

**For technical support, please call 800-352-1512 (US Customers) or +1 (858) 673-3700 (International Customers) and ask to speak with a Restech technical support specialist.**